



Open-Motion Research

User Manual

Product	Open-Motion blood-flow monitor — Research build
Application version	1.5.3
Compatible firmware	Console 1.8.1 · Sensor modules 1.8.2
Document revision	September 2026 — draft
Audience	Research and study operators

DRAFT — not a controlled document. This manual is a documentation preview written against Open-Motion Research 1.5.3 and illustrated with screenshots of that version running on bench hardware. Open-Motion Research is intended for investigational use only.

Contents

1. About the system
2. Installing, updating and starting
3. The main screen
4. Scan Settings
5. Checking contact quality
6. Running a scan
7. Reading and navigating the plots
8. Session Notes
9. Scan History
10. Settings
11. Data storage
12. Messages and errors
13. Getting help
14. Known issues in 1.5.3

1. About the system

Open-Motion is Openwater's optical blood-flow monitor. Two forehead **sensor modules** (eight cameras each) and a **console** containing a near-infrared laser measure laser-speckle images through the skin and compute, 40 times per second for every active camera:

- **BFI — Blood Flow Index:** a relative index of blood flow (white trace).
- **BVI — Blood Volume Index:** a relative index of blood volume (blue trace).
- **Mean and Contrast:** the raw speckle-image statistics that BFI and BVI are derived from, available as an alternative display.

Open-Motion Research is the research build of the Open-Motion application. It shows an orange **BETA** badge in the header and gives full control over acquisition: you choose which cameras to record, whether a scan is timed or continuous, and how the data is displayed, with one live plot per camera.

Laser safety is enforced by a hardware interlock inside the console, independently of this software.

2. Installing, updating and starting

2.1 Installing

- Download **Open-Motion-Research-Setup-1.5.3.exe** from the Open-Motion releases page on GitHub (github.com/OpenwaterHealth/openmotion-bloodflow-app/releases). The installer and the application are code-signed by Openwater.
- The installer installs the USB driver for the sensor modules, then the application, and offers a desktop shortcut. It needs administrator rights.
- **Each Windows account has its own scan history, settings and logs** (§11).
- **Upgrading from 1.5.2 or earlier:** the application starts with an empty scan history and default display preferences. Earlier scans are not deleted — they remain in `C:\ProgramData\Openwater\data` — but this version does not show them.
- A macOS build ([Open-Motion-1.5.3-macOS.dmg](#)) is published on the same page. This manual describes the Windows build.

2.2 Updating the application

About 3 seconds after launch, the application checks GitHub for a newer release (this needs internet access). If one exists, a banner appears under the header: "A new version is available: X.Y.Z" with an **Update** button and a  to hide it. The same status appears in Settings → About (§10).

Pressing **Update** downloads the new installer ("Downloading update..."), checks that it carries a valid Openwater signature, then closes the application, installs the update and starts the new version ("Installing update..."). Windows may ask for permission to install. An installer without a valid Openwater signature is refused and nothing is installed. You can always install a release manually instead (§2.1).

2.3 Updating device firmware

When newer firmware is available for a connected device, a banner reads "Device firmware update available". **View** opens Settings, where the device's row in the About card shows an **Update** button. Pressing it asks for confirmation:

"Console FW will reboot into DFU mode and be re-flashed. Do not unplug it until this completes."

Progress is shown under the firmware rows. Do not unplug or power off the device while it updates, and do not update firmware while a study is in progress.

2.4 Starting the application

Launch **Open-Motion Research** from the Start menu or the desktop shortcut. The application is a single signed file that unpacks itself on launch, so the window takes a few seconds to appear.

- Only one copy of Open-Motion can run at a time. A second launch shows *"Another instance of the application is already running. Please close the existing instance before opening a new one."*
- The application connects to the console and the sensor modules **automatically** whenever they are plugged in and powered on. There is no Connect button.
- About 12 seconds after launch the application checks what it found. If a device is missing, a yellow message appears in the bottom-right corner:
 - *"Console not detected. Check the console USB cable and power, then reconnect."*
 - *"Sensor not detected. Check the sensor USB cable and power, then reconnect."*
 - *"System not found. Check that the console and sensor are connected and powered on."* (nothing was found)
- **Sample dataset:** if no device at all is connected at that point, a dialog titled *"No console detected"* offers to open a sample recording — *"You can explore real BFI/BVI traces without any hardware connected."* Choose **Open sample dataset** to load it into the plot viewer, or **Not now**. The offer is made once per launch.
- After a sensor module connects, the system initializes it for a few seconds. The **Start** label stays dimmed and cannot be pressed until initialization finishes.

2.5 Connection states


The round badge at the top of the left toolbar shows the system state:

Badge	Label	Meaning
Grey circle, chain-link icon	Disconnected	The console or both sensor modules are not connected. Scanning is disabled.
Green circle, dimmed label	Start	Connected; a sensor module is still initializing. Wait a few seconds.
Green circle, play icon	Start	Ready to scan.
Yellow circle	Start	Start was pressed; the system is finishing the previous step before the scan begins (normally a few seconds).
Red circle, stop icon	Stop	A scan is running. Press to stop it.

2.6 Console indicator light

Console light	Meaning
Solid green	Powered and idle.
Solid blue	The laser is active (during a check or a scan).
Blinking blue	An error. The application shows an error dialog (§12.3); the light returns to green once the dialog is dismissed. If the light blinks with no dialog in the application, power-cycle the console.

3. The main screen

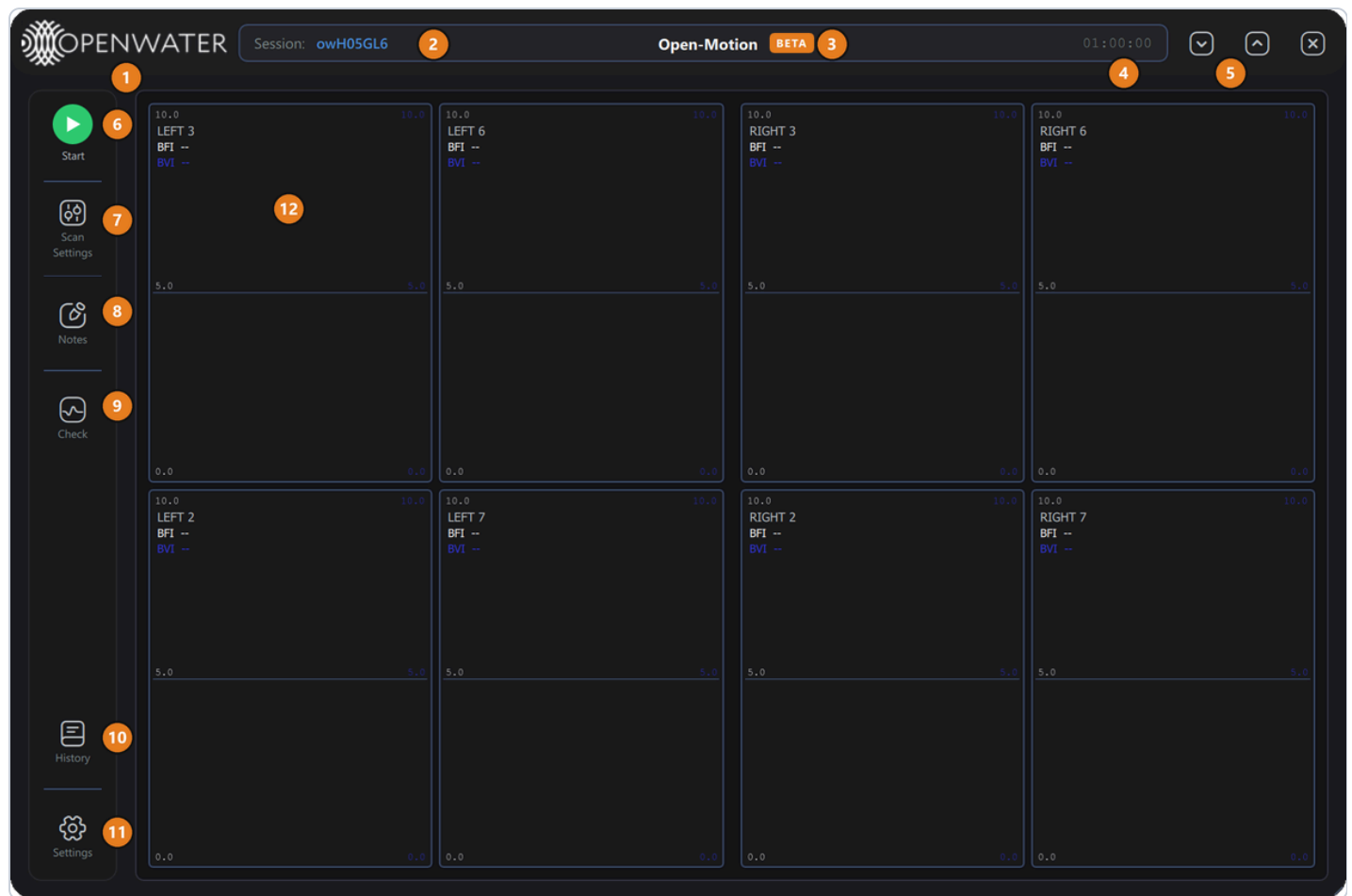

#	Element	What it does
1	Openwater logo	Branding.
2	Session	The session label (e.g. <code>owH05GL6</code>). A new label is generated at each launch; set your own in Scan Settings (§4). Stored with every scan.
3	BETA badge	Identifies the Research build.
4	Scan clock	Idle: the configured scan length (e.g. <code>01:00:00</code>) or <code>Continuous</code> . Scanning: <code>elapsed / total</code> in green, or <code>Continuous</code> and the elapsed time.
5	Window controls	Minimize, maximize/restore, and close.
6	Start / Stop badge	Starts the scan immediately, or stops the running scan (§2.5).
7	Scan Settings	Session label, camera selection and scan length (§4). Disabled during a scan.


#	Element	What it does
8	Notes	Opens Session Notes (§8). Always available.
9	Check	Runs a contact-quality check (§5). Disabled during a scan.
10	History	Opens Scan History — review, replay, export or delete past scans (§9). Disabled during a scan.
11	Settings	Opens Settings (§10). Disabled during a scan.
12	Plot grid	One plot per active camera, labeled <code>LEFT n</code> / <code>RIGHT n</code> (§7).

Window behavior:

- **Move** the window by dragging the header bar.
- **Resize** with the diagonal grip in the bottom-right corner (minimum 800 × 600).
- **Close while busy:** if you press `×` during a scan or check, the application does not exit immediately. A yellow message such as *"Scan in progress. Click X again to cancel and exit."* appears. Press `×` again within 5 seconds to cancel the work and exit.

4. Scan Settings

Press **Scan Settings** to set up the next scan. Closing the dialog (**X** , **Esc** , or clicking outside it) applies the settings.

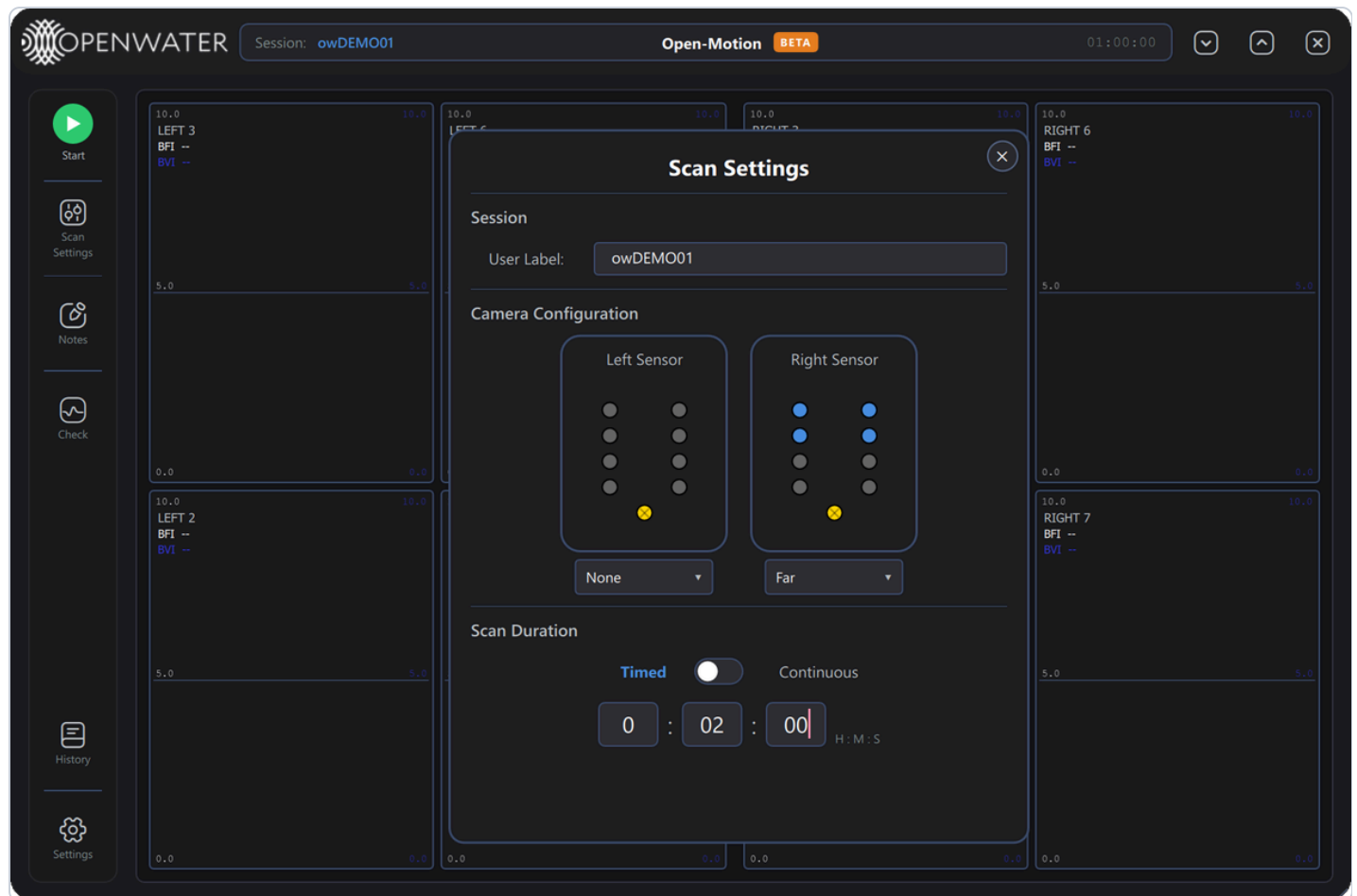

*Scan Settings with the label **owDEMO01** , no cameras on the left sensor, the Far pattern on the right sensor, and a 2-minute timed scan.*

4.1 Session


User Label names the session. Letters and digits are kept, converted to upper case and prefixed with **ow** — typing **demo01** gives **owDEMO01** . The label appears in the header, in Scan History, and in every scan's full label (e.g. **20260925_222432_owDEMO01** , with the scan's start time in UTC). The label stays in effect until you change it or restart the application.

Use a study or subject code, never a name or other patient-identifiable information.

4.2 Camera Configuration

Each sensor diagram shows the module's eight cameras: cameras 1 to 4 down the left column, 8 to 5 down the right column, and the laser aperture (yellow ⊗) at the bottom. **Blue** cameras will be recorded; **grey** cameras will not. Hover over a camera to see its number. A disconnected sensor's diagram is dimmed and cannot be changed.


Choose a pattern for each sensor from the list under its diagram:

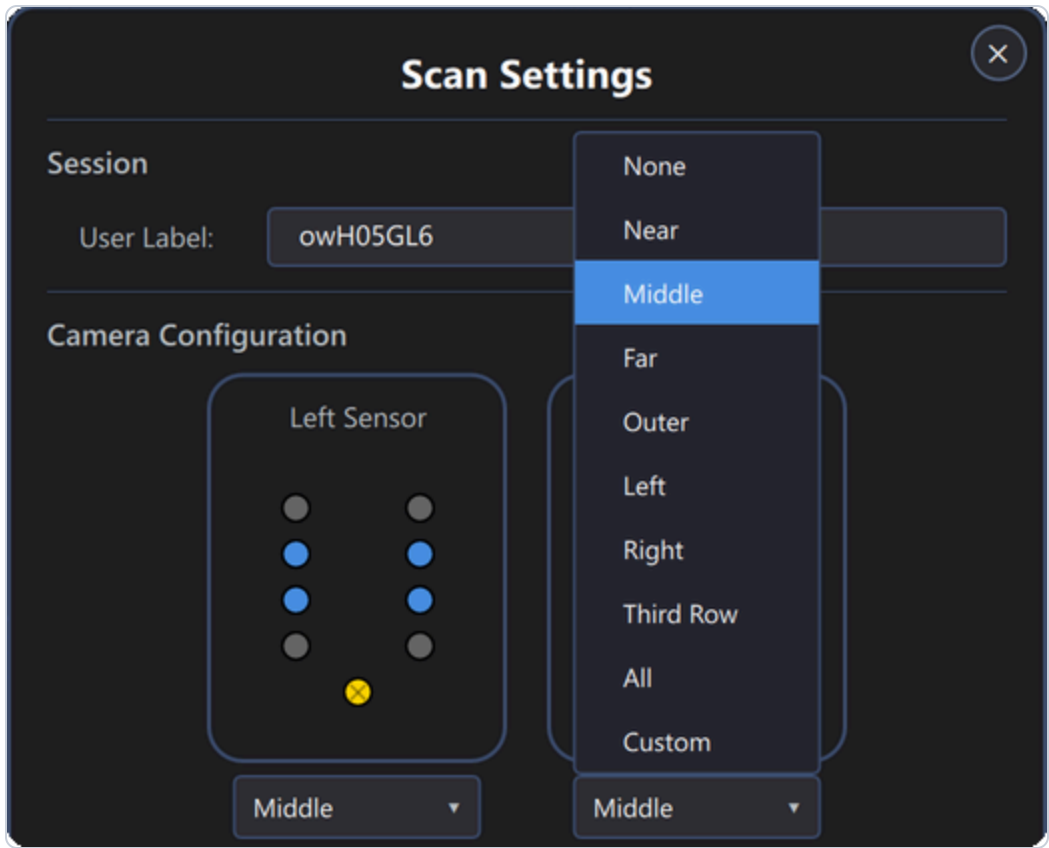

Pattern	Cameras recorded
None	none — the sensor is not used
Near	2, 4, 5, 7
Middle	2, 3, 6, 7 (default)
Far	1, 2, 7, 8
Outer	1, 4, 5, 8
Left	1, 2, 3, 4
Right	5, 6, 7, 8
Third Row	2, 7
All	all eight

Pattern	Cameras recorded
Custom	any combination — see below

Custom: a ring appears around each camera and clicking a camera turns it on or off. Custom starts from the cameras currently shown.

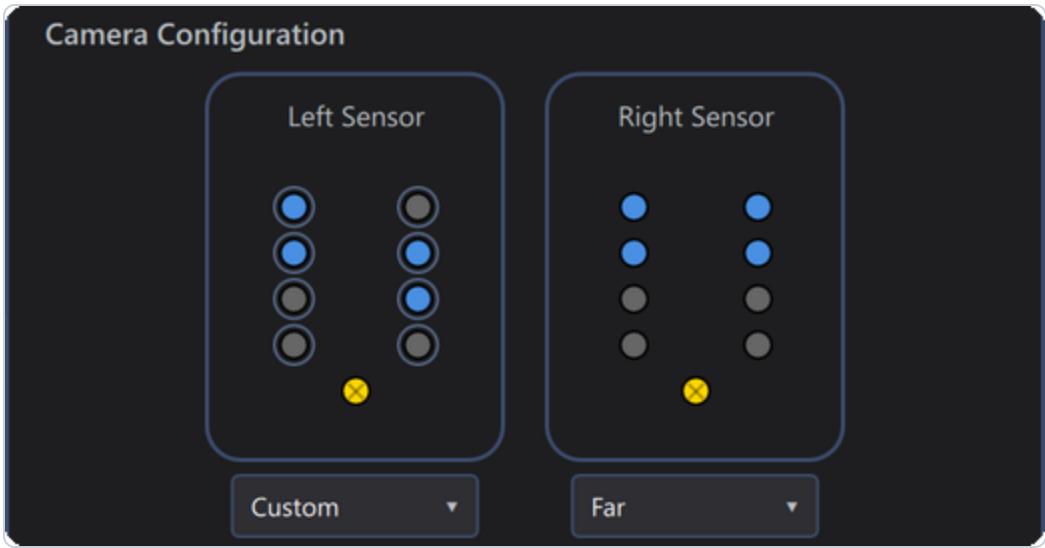


Custom mode on the left sensor: each camera can be clicked to include or exclude it.

More cameras give more coverage but denser plots and larger recordings. The selection is kept if a sensor briefly disconnects and reconnects. The patterns selected when the application starts come from Settings → Default Camera Configuration (\$10).


4.3 Scan Duration

Control	What it does
Timed / Continuous switch	Timed: the scan stops by itself after the set time. Continuous: the scan runs until you press Stop ("Scan will run indefinitely until stopped."), up to a maximum of 12 hours.
H : M : S	The timed length (default 1:00:00, up to 99:59:59). A length of 0:00:00 is not allowed; it is reset to 1 minute with the message "Scan duration cannot be 0 seconds — reset to 1 minute."

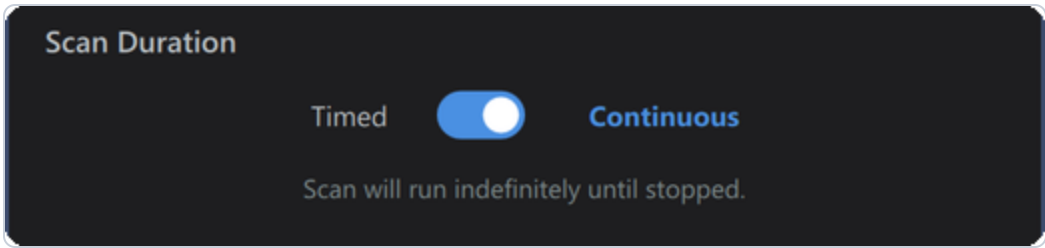

5. Checking contact quality

Press **Check** (not available during a scan) to check sensor contact without recording. The laser is on during the check, which takes about 10 seconds. The check tests **all eight cameras** on each connected sensor, whatever is selected in Scan Settings.

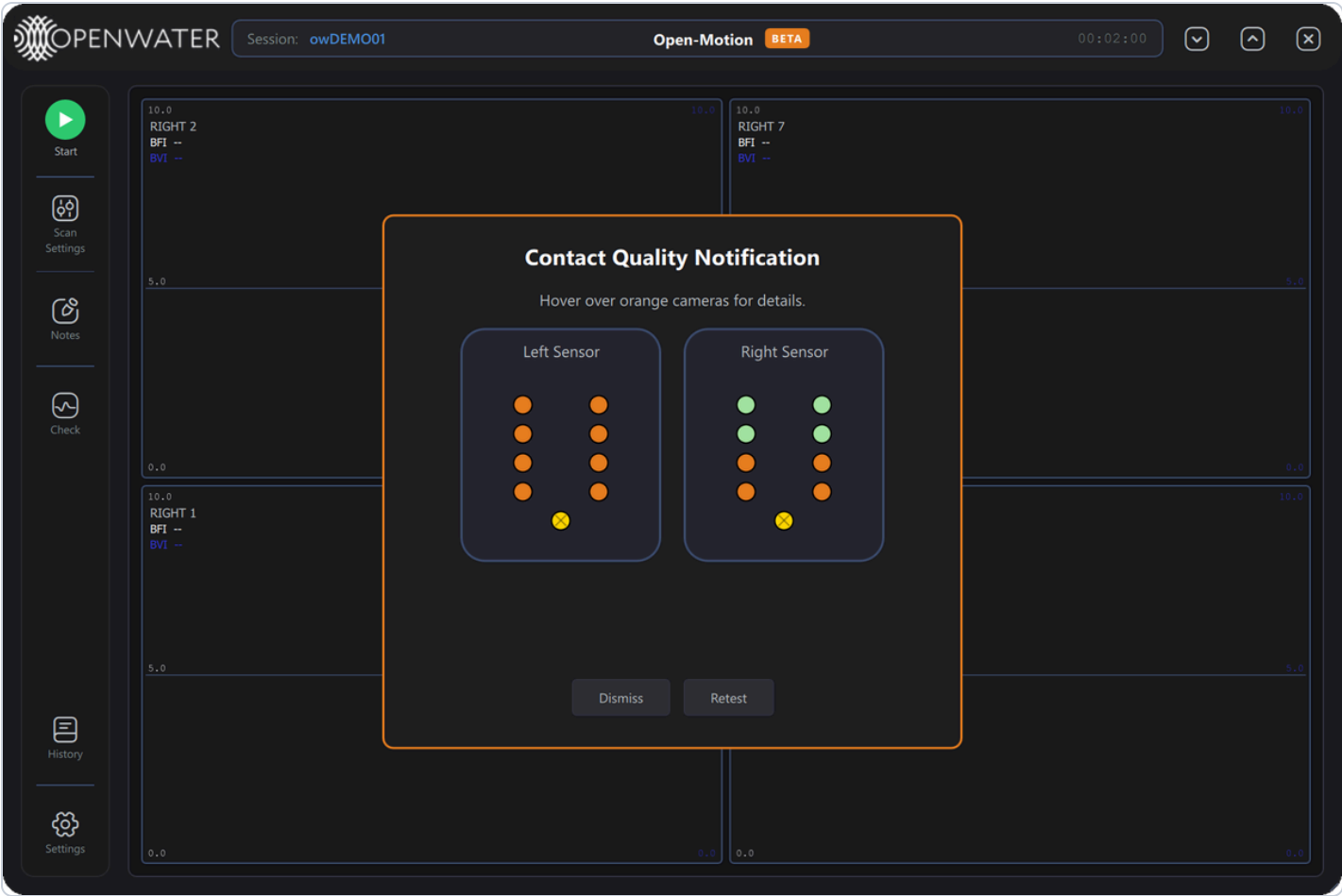

A check with problems: the left sensor was not in contact, and the right sensor had good contact on its top four cameras only.

Border	Title	Meaning
Green	Good signal quality	All cameras report acceptable ambient light and contact levels.
Orange	Contact Quality Notification	One or more cameras have a problem. Hover over an orange camera for the reason.
Red	Contact check failed	The check could not be completed; the reason is shown in red.


Camera dots: **green** = good contact · **orange** = problem · **grey** = not evaluated. Hovering over a camera shows its label (L = left, R = right, then the camera number) and the problem:

Reason	What to do
Poor sensor contact	The camera is not receiving enough light back from the skin. Re-seat the module flat against the skin and check for hair or debris.
Ambient light detected	Room light is reaching the camera. Improve the seal against the skin or shade the sensors.

Dismiss closes the dialog; **Retest** runs the check again after you adjust the sensors. A scan does not run the check for you — run it before pressing Start when contact matters.

6. Running a scan


1. Set up the scan in **Scan Settings** (§4). The defaults (Middle pattern, 1-hour timed scan) are a reasonable starting point.
2. Run **Check** (§5) and re-seat the sensors until the cameras you plan to record are green.
3. Press **Start**. The scan starts immediately: the badge may turn yellow for a moment, then red (**Stop**), and the header clock starts counting.

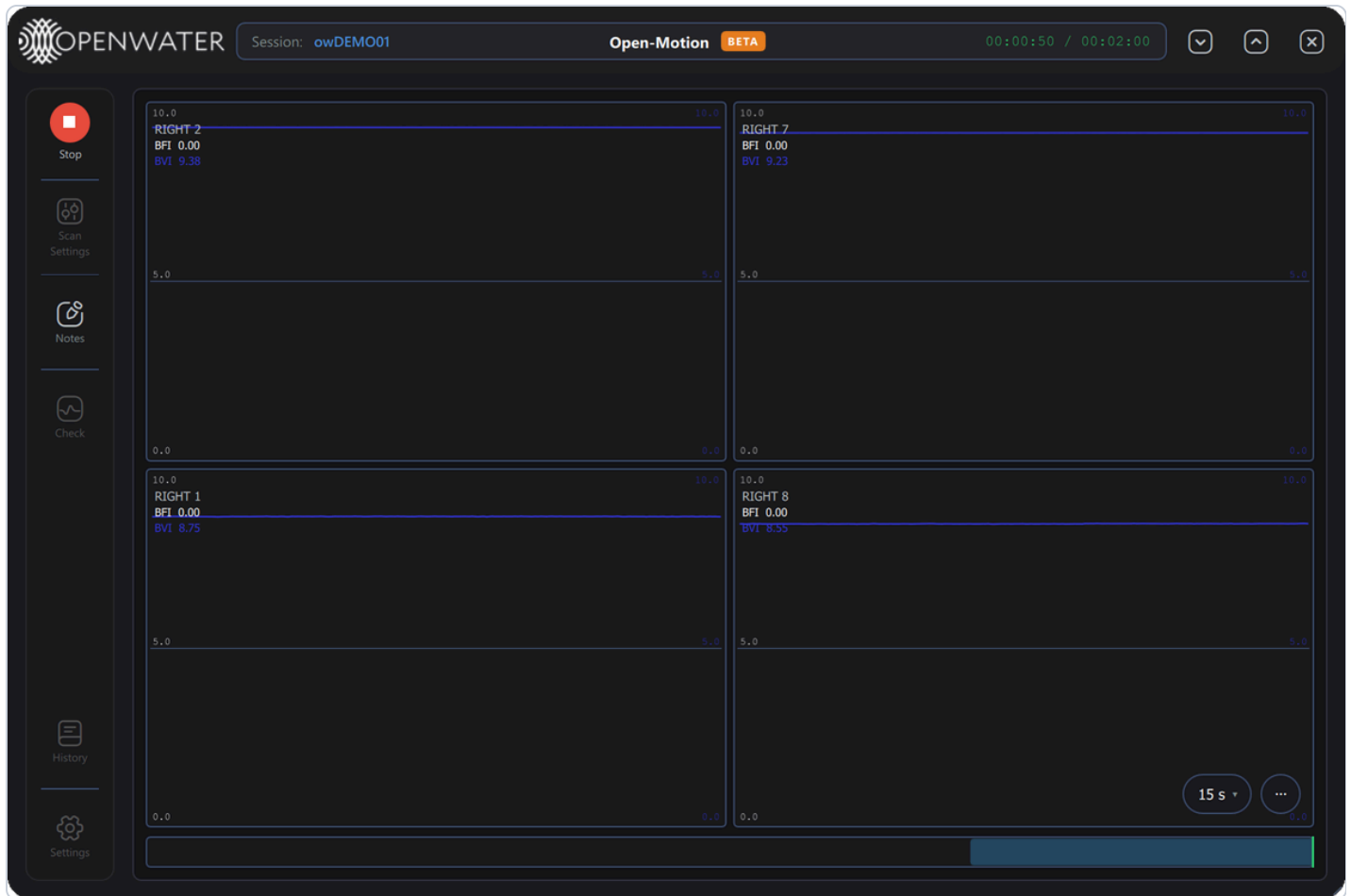

A 2-minute timed scan on four cameras of the right sensor, 26 seconds in. Bench capture — the sensors were not on a subject, so the values are not physiological.

During the scan:

- Each plot shows one camera's BFI (white) and BVI (blue) traces, with the latest values under its label (§7).
- **Scan Settings**, **Check**, **History** and **Settings** are disabled. **Notes** stays available.
- Press **Space** to open Session Notes with a timestamp already inserted:

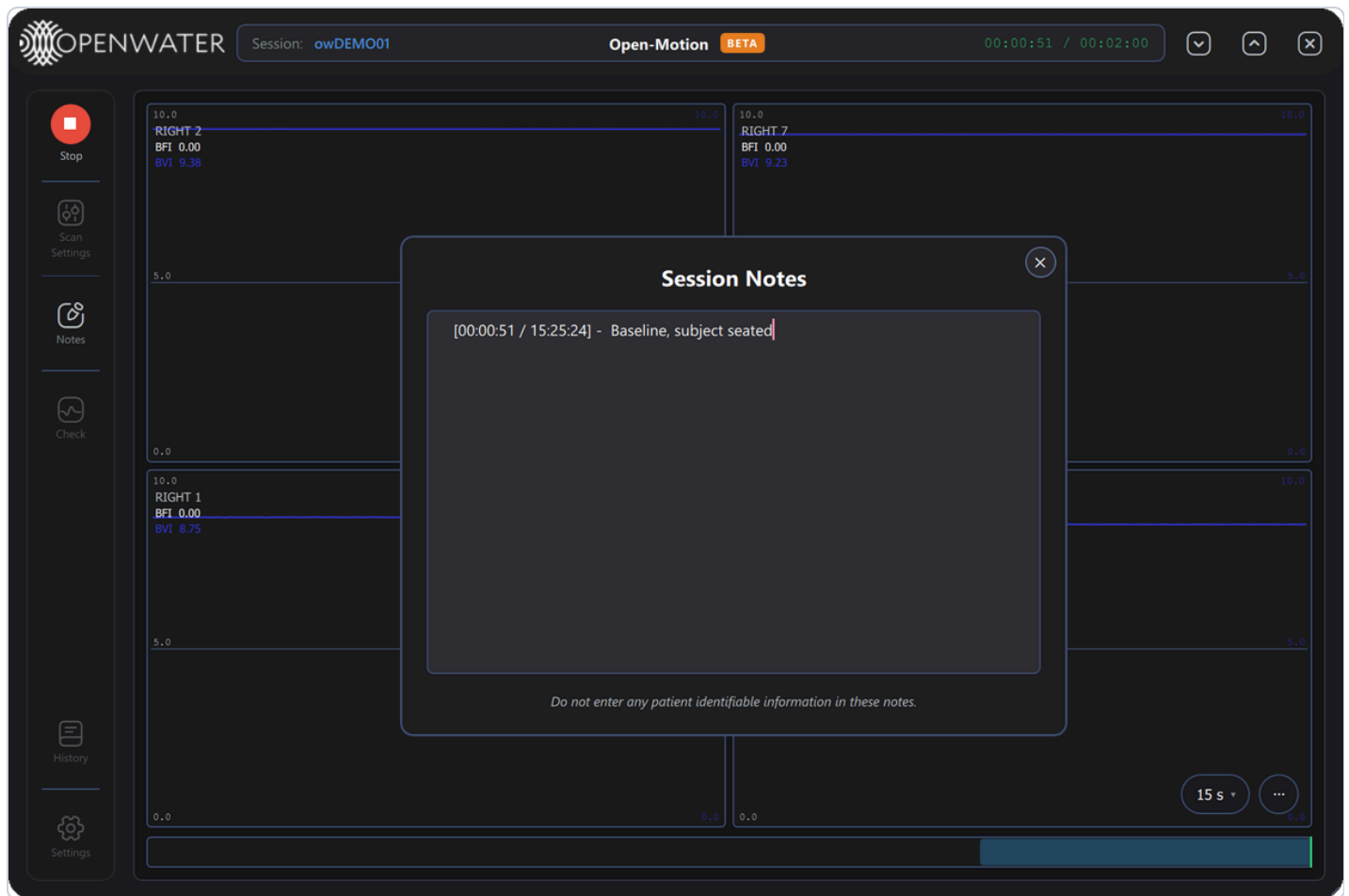

Space during a scan inserts `[elapsed / clock time] -` — type the observation and close the window with `X` to save it.

Contact warnings during a scan. If contact degrades on a recorded camera, the contact-quality dialog opens on top of the plots. Only the recorded cameras are evaluated; the others are grey. The scan keeps recording while the dialog is shown.

Button	Action
Stop scan	End the scan now.
Continue	Return to the plots. Enabled only after every problem has stayed clear for a few seconds — re-seat the sensor and wait for the border to turn green and the message "All contact quality issues are currently inactive. You may dismiss."


Camera messages. Short messages in the bottom-right corner report individual cameras:

- "Camera RIGHT 2 connection lost at 00:01:23 — last temp 41.2°C" — that camera has delivered no data for more than 2 seconds, and its plot shows **CONNECTION LOST**. The scan continues on the remaining cameras. Check the sensor module's cable and power.
- "Camera RIGHT 2 reconnected at 00:01:40" — data from that camera resumed.

- "Camera RIGHT 2 temperature ... — above 110°C threshold. Check the airflow around the sensor module." — a camera is running hot. Make sure nothing covers the module's vents.

If every camera stops delivering data, the scan is stopped automatically (E-303, §12.4).

Ending the scan. A timed scan stops by itself; otherwise press the red **Stop** badge (or **Stop scan** in a contact warning). After a few seconds, **Session Notes** opens with a closing line — "Scan completed — duration: 00:02:00" for a timed scan that ran to the end, "Scan stopped — duration: ..." when stopped early.

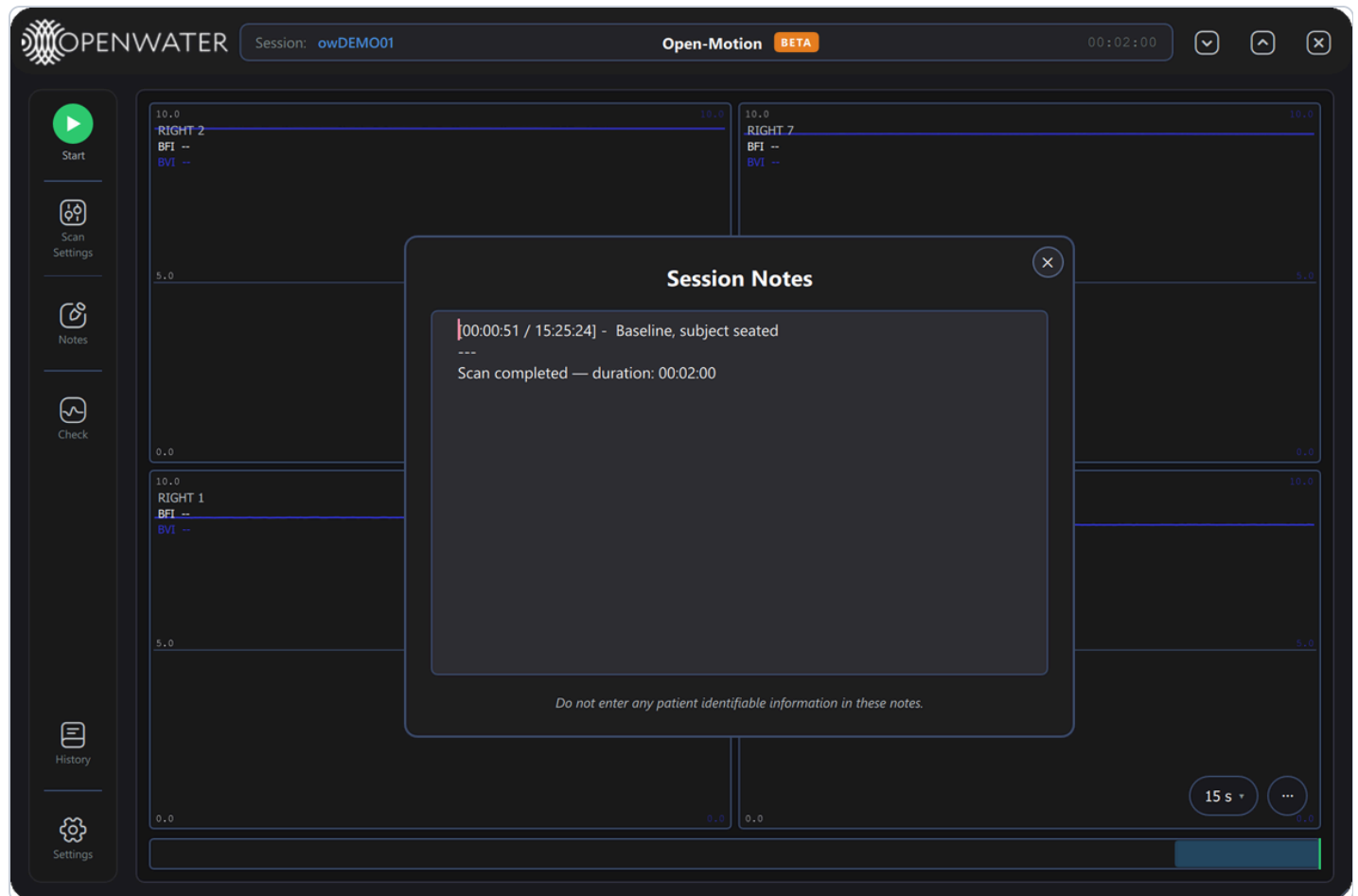

- If gaps longer than 1 second were detected in the data, a "Data gaps (> 1.0s): ..." line is added.
- All scan data is saved automatically to the scan database (§11); raw CSV files are also written if enabled in Settings (§10).
- If a scan ends unexpectedly, a message reports what was kept:
 - "Scan ended unexpectedly — partial data was saved. The final segment could not be dark-corrected and was discarded." (yellow)
 - "Scan ended unexpectedly and no data was recorded (the device may have disconnected mid-scan). This scan was not saved." (red; stays on screen until dismissed)

7. Reading and navigating the plots

7.1 The plot grid

There is one plot per recorded camera. The grid mirrors the sensor modules: the left module's cameras on the left, the right module's on the right, and each row pairs the cameras that sit opposite each other — 4 and 5, 3 and 6, 2 and 7, 1 and 8 from top to bottom. Rows with no recorded camera are left out.


Each plot shows:

- the camera label (`LEFT 3` , `RIGHT 7` , ...) and its latest values under it;
- the traces: BFI (white) and BVI (blue), or Mean (green) and Contrast (purple);
- axis labels: the left scale belongs to the first trace (BFI or Mean), the right scale to the second (BVI or Contrast);
- **CONNECTION LOST** while that camera is not delivering data during a live scan.

BFI and BVI values are limited to a 0–10 display range: a value below 0 shows as `0.00` and a value above 10 as `10.00`. `--` means there is no valid reading.

7.2 Display options — the `...` menu

The `...` button at the bottom-right of the plots opens the display options:

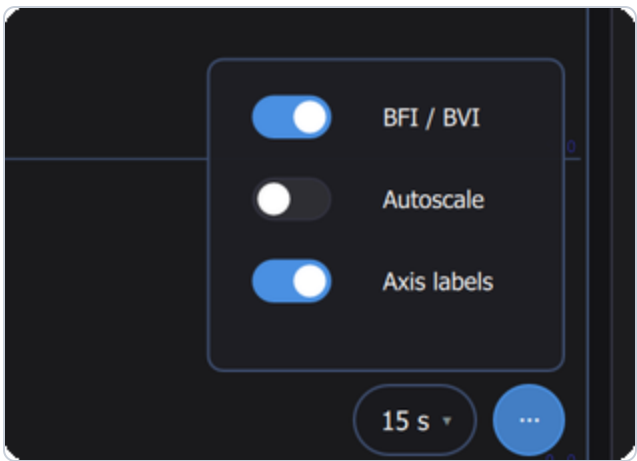


Switch	What it does
BFI / BVI ↔ Mean / Contrast	Switches every plot between BFI/BVI and the raw Mean/Contrast statistics.
Autoscale	Fits the vertical axis to the data instead of using the Manual Plot Bounds from Settings.
Global scale ↔ Per-plot scale	Shown while Autoscale is on. Global: all plots share one range. Per-plot: each plot fits its own data.

Switch	What it does
Axis labels	Shows or hides the numbers on the vertical axes.

These switches are the same settings as in Settings → Realtime Plot Display and stay set across restarts.



Autoscale on: each axis fits the data, so small changes become visible.

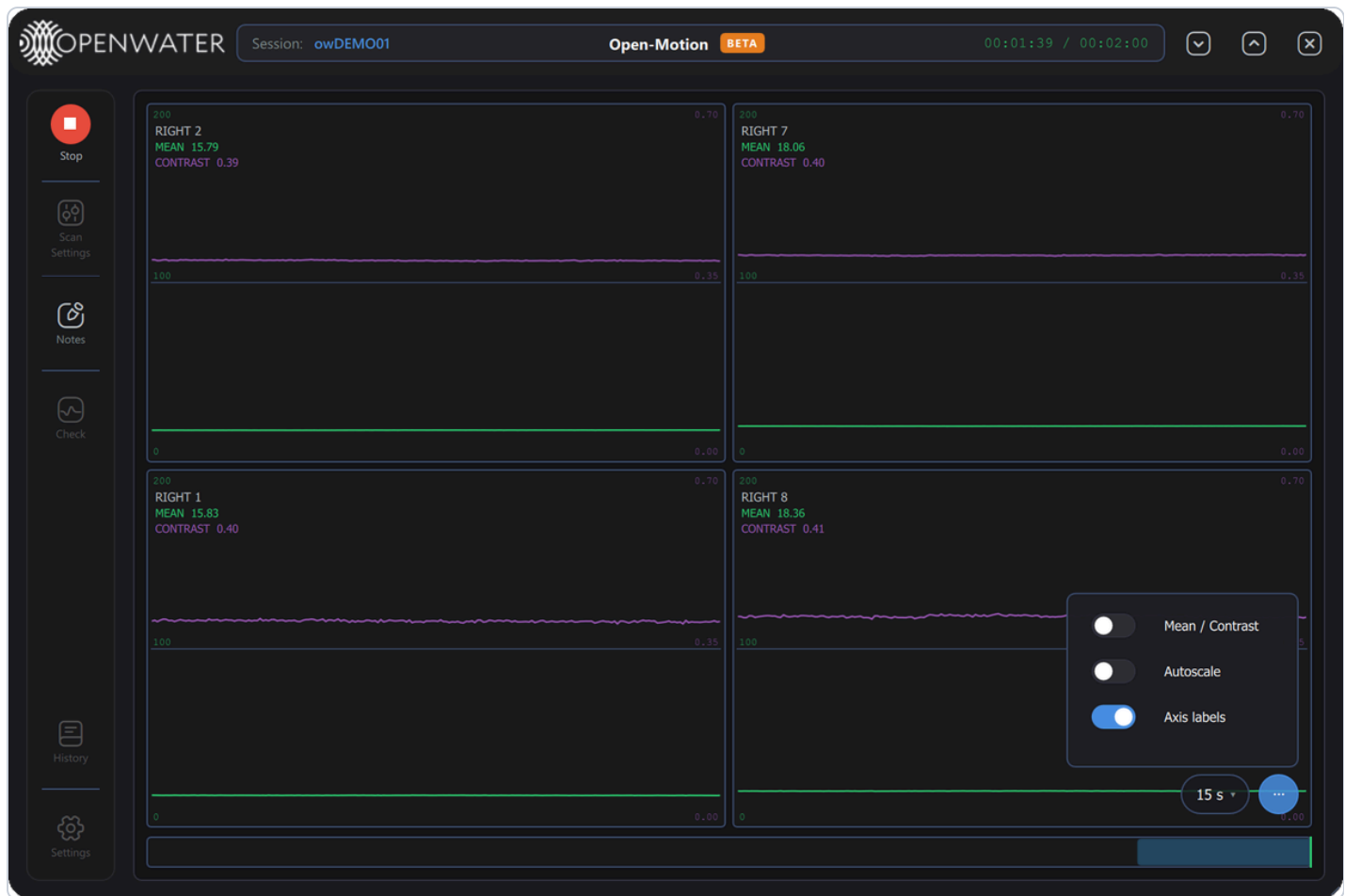

Mean / Contrast display: Mean (green) and speckle Contrast (purple) for each camera.

7.3 Navigation

The plot area works like a video recorder: you can rewind during or after a scan without losing the live recording.

Control	Where	What it does
Window-length pill (e.g. 15 s ▾)	Bottom-right of the plots	Chooses how much time is visible: 5 s, 15 s, 30 s, 1 min or 5 min. The plots open at 15 s.
Timeline bar	Under the plots	Shows the whole scan. Drag the highlighted window to move through time, or click the bar to jump. Blue = following the live edge; orange = paused in the past.
Drag on a plot	Plot area	Pans back or forward in time (stops following live).
Mouse wheel on a plot	Plot area	Zooms the time window in or out around the pointer.
Hover	Plot area	Shows a readout box (top-right) with the time and every camera's values under the pointer.
• Back to live	Top-right of the plots	Appears when you have moved into the past during a scan. Click to return to the live edge.

Keyboard shortcuts (active when the plot area has keyboard focus — click a plot first):

Key	Action
← / →	Pan 1 second back / forward (hold Shift to pan a full window)
+ or ↑	Zoom in
- or ↓	Zoom out
0	Reset the zoom and return to live
Home / End	Jump to the start of the scan / return to the live edge
Esc	Release keyboard focus from the plots
Space	Open Session Notes with a timestamp (during a scan)


8. Session Notes

- Open the notes with the **Notes** button at any time. They open automatically after every scan, and **Space** during a scan opens them with a timestamp line ([elapsed / clock time] -) ready for you to type after.
- Notes are attached to a scan. During a scan they belong to that scan; after a scan they edit the scan just recorded. Each new scan starts with empty notes, so text typed before a session's first scan is cleared when that scan starts.
- Notes are saved when the window closes — with `×`, `Esc`, or by clicking outside the window. "Note saved." confirms it. **Always close the Notes window after typing** so the text is saved.
- Saved notes appear, read-only, in Scan History (§9). They are kept in the scan database and are not included in exported CSV files.

Do not enter patient-identifiable information in the notes. The same reminder is shown at the bottom of the Notes window.

9. Scan History

Press **History** (available when no scan is running) to review past scans.

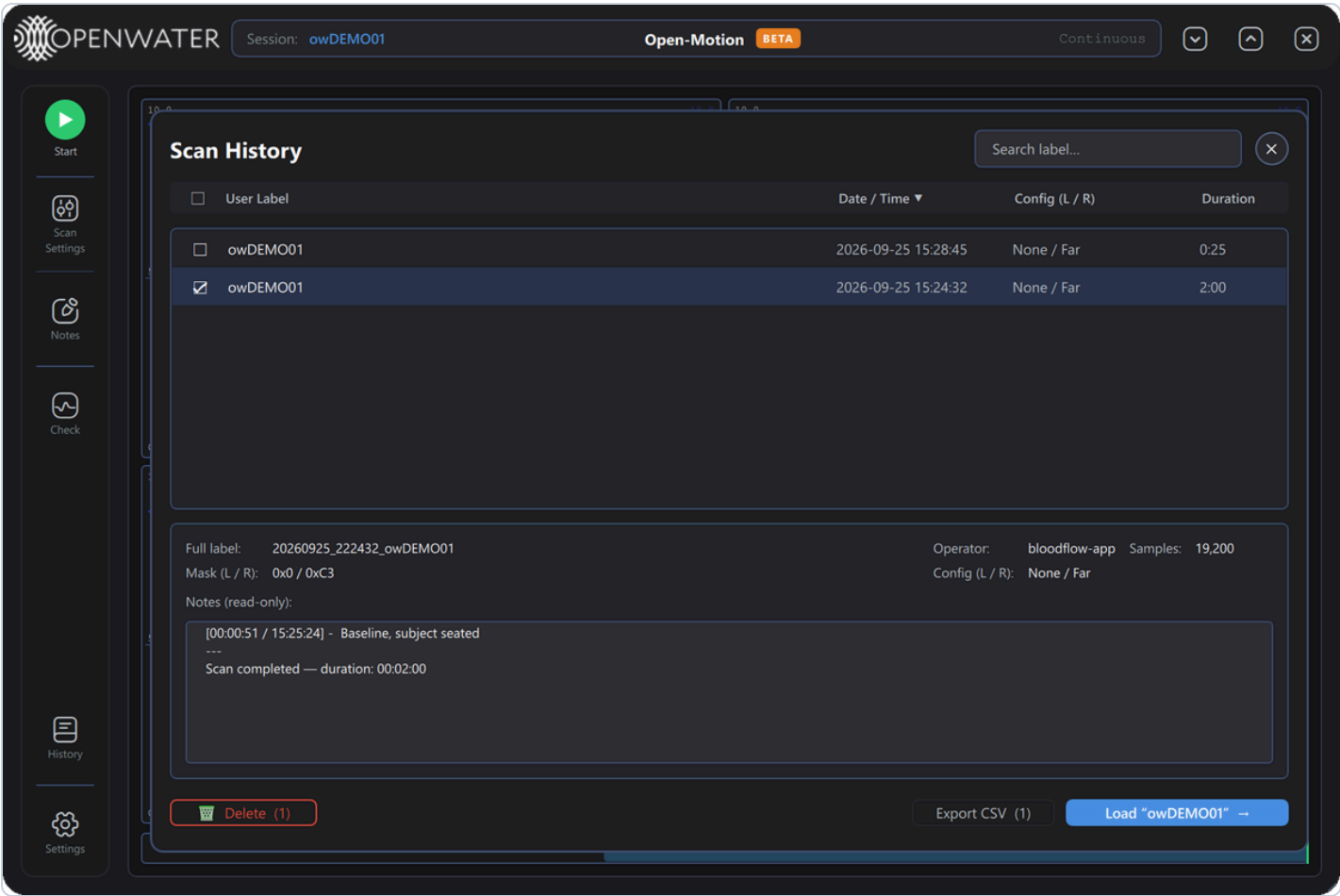

Column / control	What it shows or does
User Label	The session label the scan was recorded under. Scans from one session share a label; use the date and time to tell them apart.
Date / Time	When the scan started, in the computer's local time.
Config (L / R)	The camera pattern used on each side (<code>None</code> , <code>Middle</code> , <code>Far</code> , ...). A Custom selection shows as its camera mask in hexadecimal, e.g. <code>0x5B</code> .
Duration	Scan length in minutes:seconds. <code>—</code> means the scan was interrupted and has no recorded end.
Search label...	Filters the list by label.
Column headers	Click to sort; click again to reverse (▲ / ▼).

Column / control	What it shows or does
Row click	Selects the scan and shows its details: full label, operator, number of samples, camera masks (hexadecimal, left / right), configuration and notes (read-only).
Checkboxes	Mark several scans for Delete or Export. The header checkbox selects all.
Export CSV (N)	Exports the checked scans to CSV (§9.1).
Load "label" →	Opens the selected scan in the plot viewer for replay (§9.3).
🗑 Delete (N)	Permanently deletes the checked scans (§9.2).
×	Closes History.

9.1 Exporting


One checked scan opens a Save dialog (default name `<full_label>_export.csv`); several ask for a destination folder and write one file per scan. Each file has one row per frame: the frame number, the time in seconds, and BFI, BVI, Mean and Contrast for every recorded camera, plus a quality flag per camera. Notes are not included. Interrupted scans cannot be exported ("*Interrupted scans can't be exported.*"). A message confirms the export ("*Exported to ...*" or "*Exported N scan(s) to ...*").

9.2 Deleting

Deleting is permanent and password-protected: pressing **Delete** opens a *Confirm Delete* prompt, and the scans are removed only after the password issued by Openwater is entered.

9.3 Replay


After **Load**, the plots show the recorded scan with a **Viewing** badge naming it (label · date and time). All plot features (§7) work the same way on a replay. When it is shown, the red ← **Back to live scan** pill returns to the live view. A scan recorded in clinical mode replays in the clinical layout (one averaged plot per side).

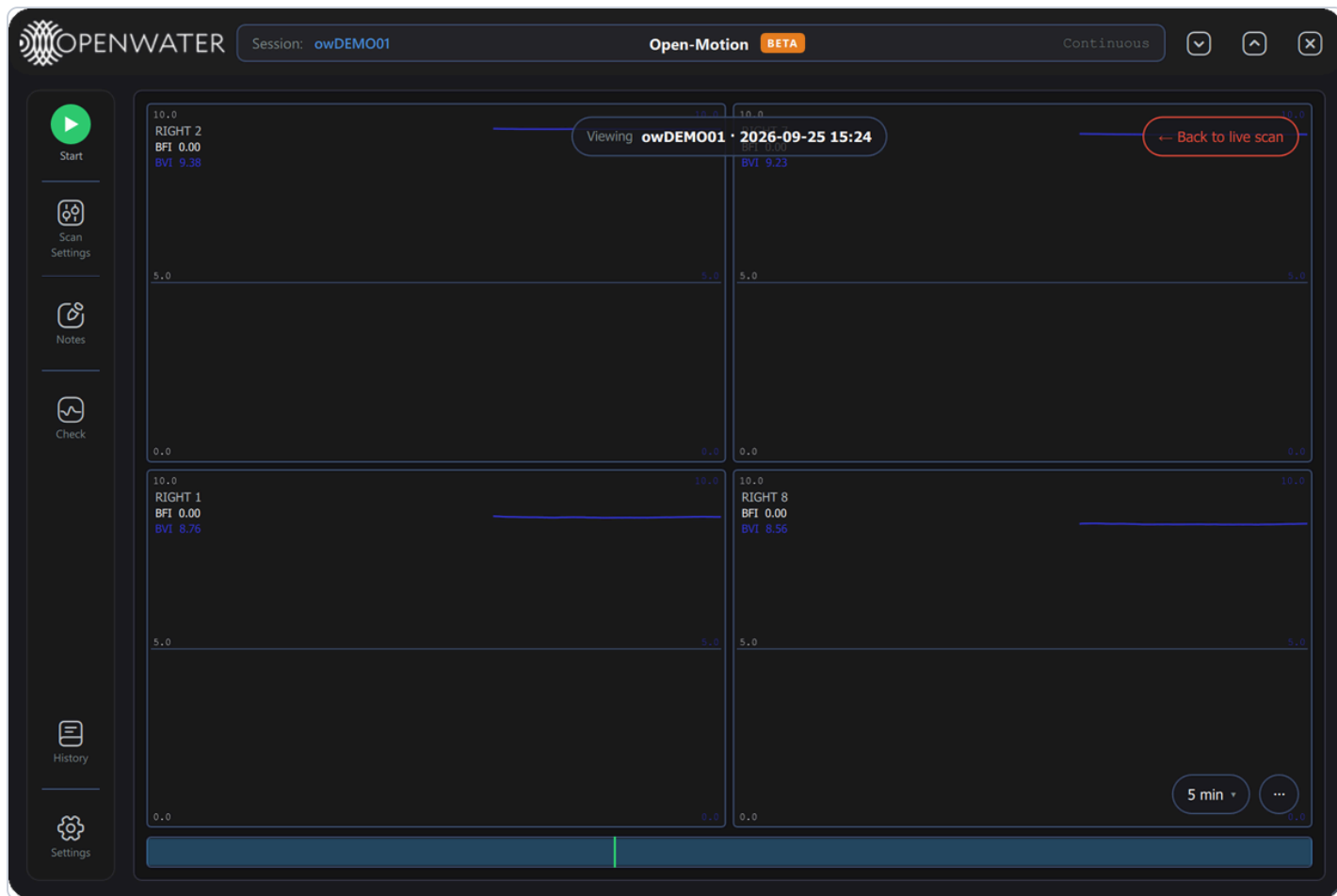

The 2-minute scan loaded from History, with the window set to 5 minutes so the whole scan fits.


10. Settings

Press **Settings** (available when no scan is running). Changes apply immediately and are saved when you close Settings with **X**, **Esc**, or by clicking outside it.

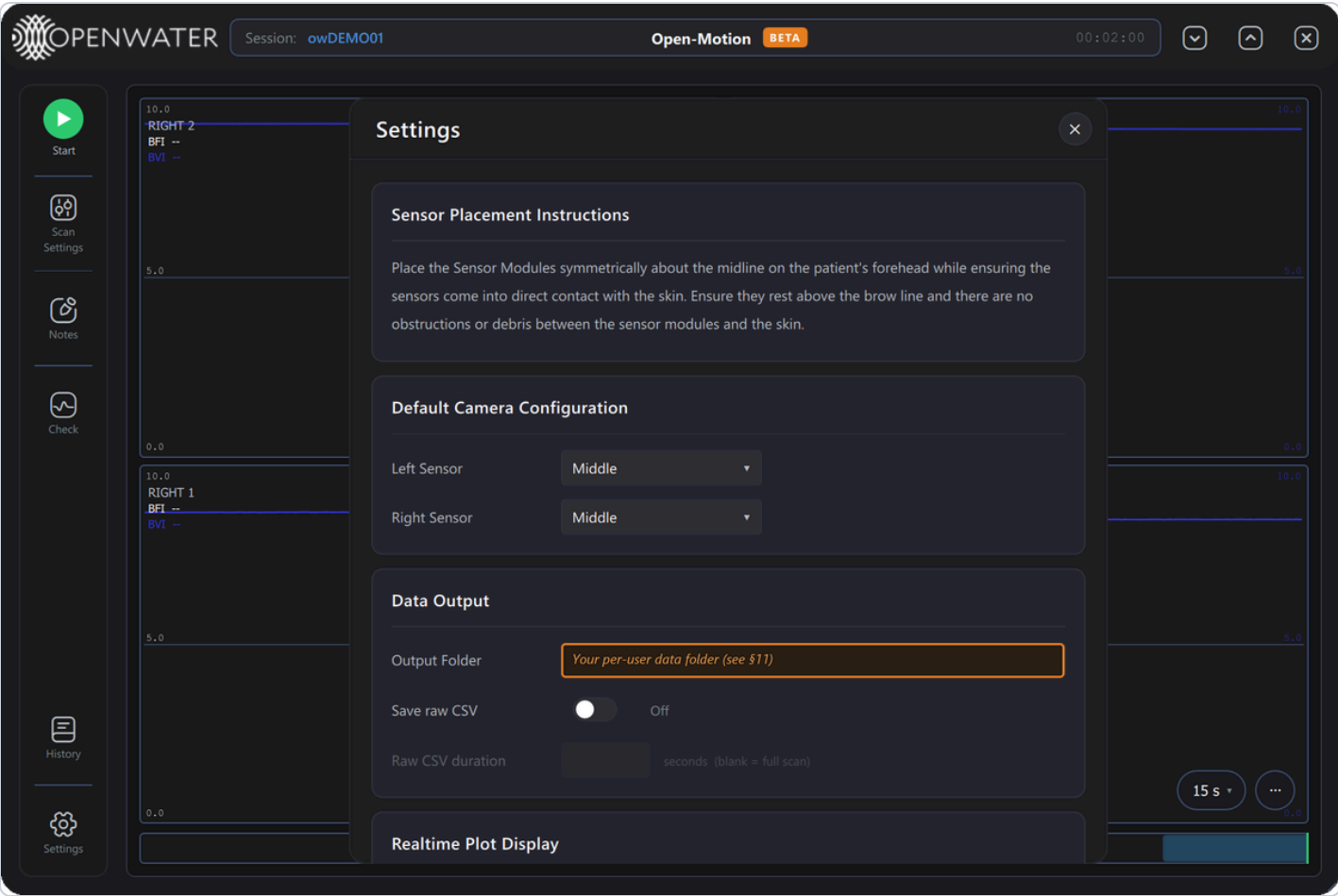

Settings, top. The Output Folder shows your data folder (§11).


Card	Contents
Sensor Placement Instructions	Place the sensor modules symmetrically about the midline on the forehead, in direct contact with the skin, above the brow line, with no obstructions or debris.
Default Camera Configuration	The Left Sensor / Right Sensor patterns (§4.2) selected when the application starts. Changes take effect at the next launch; Scan Settings changes the current session.
Data Output	<i>Output Folder</i> : this Windows account's data folder, for reference (it cannot be changed). <i>Save raw CSV</i> : also write the raw camera histograms of each scan to CSV files (§11). <i>Raw CSV duration</i> : limit raw output to the first N seconds of each scan; leave blank for the whole scan.

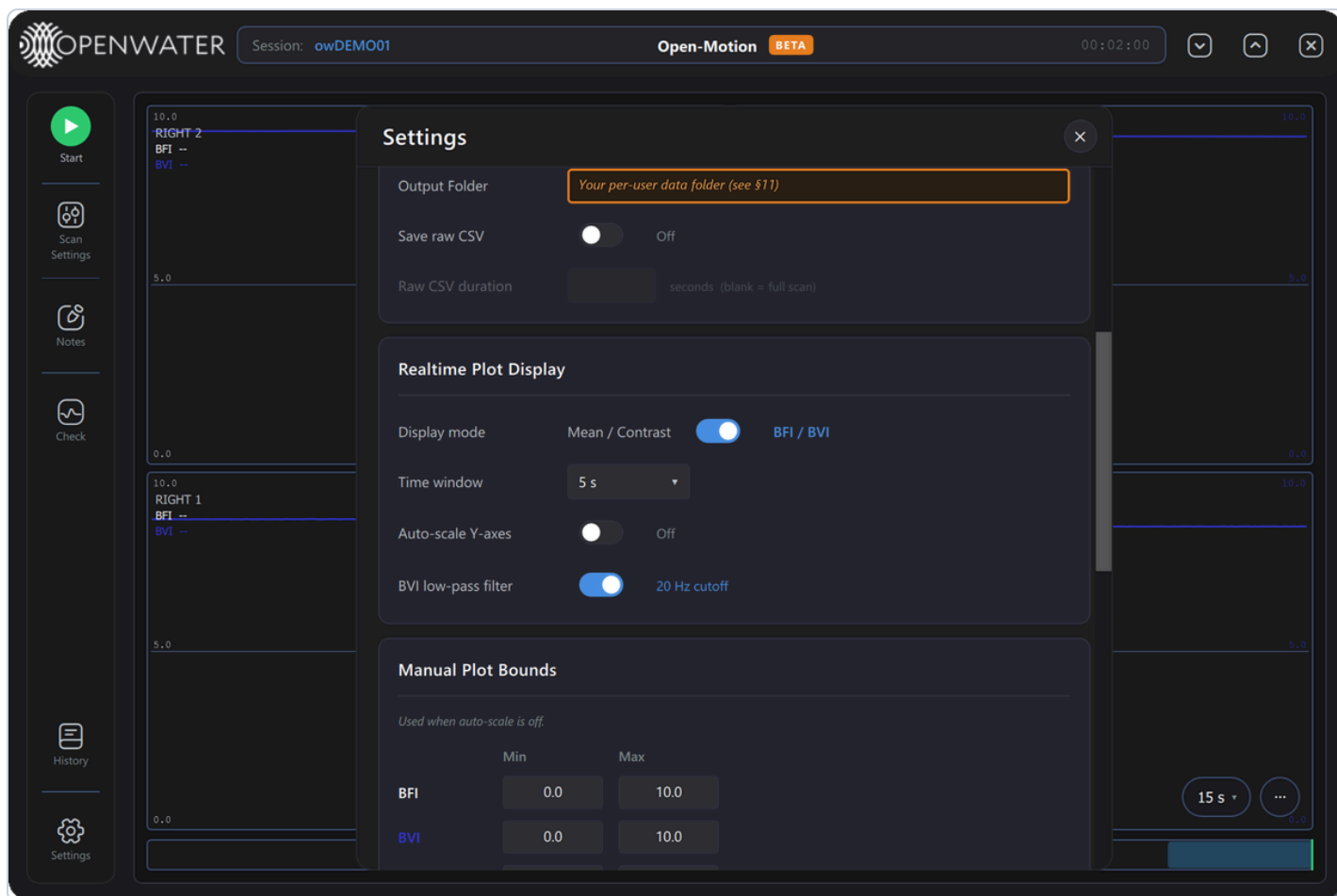


Card	Contents
------	----------

Realtime Plot Display	
------------------------------	--

	<i>Display mode</i> (Mean / Contrast ↔ BFI / BVI) and <i>Auto-scale Y-axes</i> — the same switches as the  menu (§7.2). <i>Time window</i>: 3, 5, 15 or 30 s; to change how much time the plots show, use the window-length pill (§7.3). <i>BVI low-pass filter</i>: smooths the BVI trace (20 Hz cutoff, on by default); takes effect immediately, including during a scan.
--	--

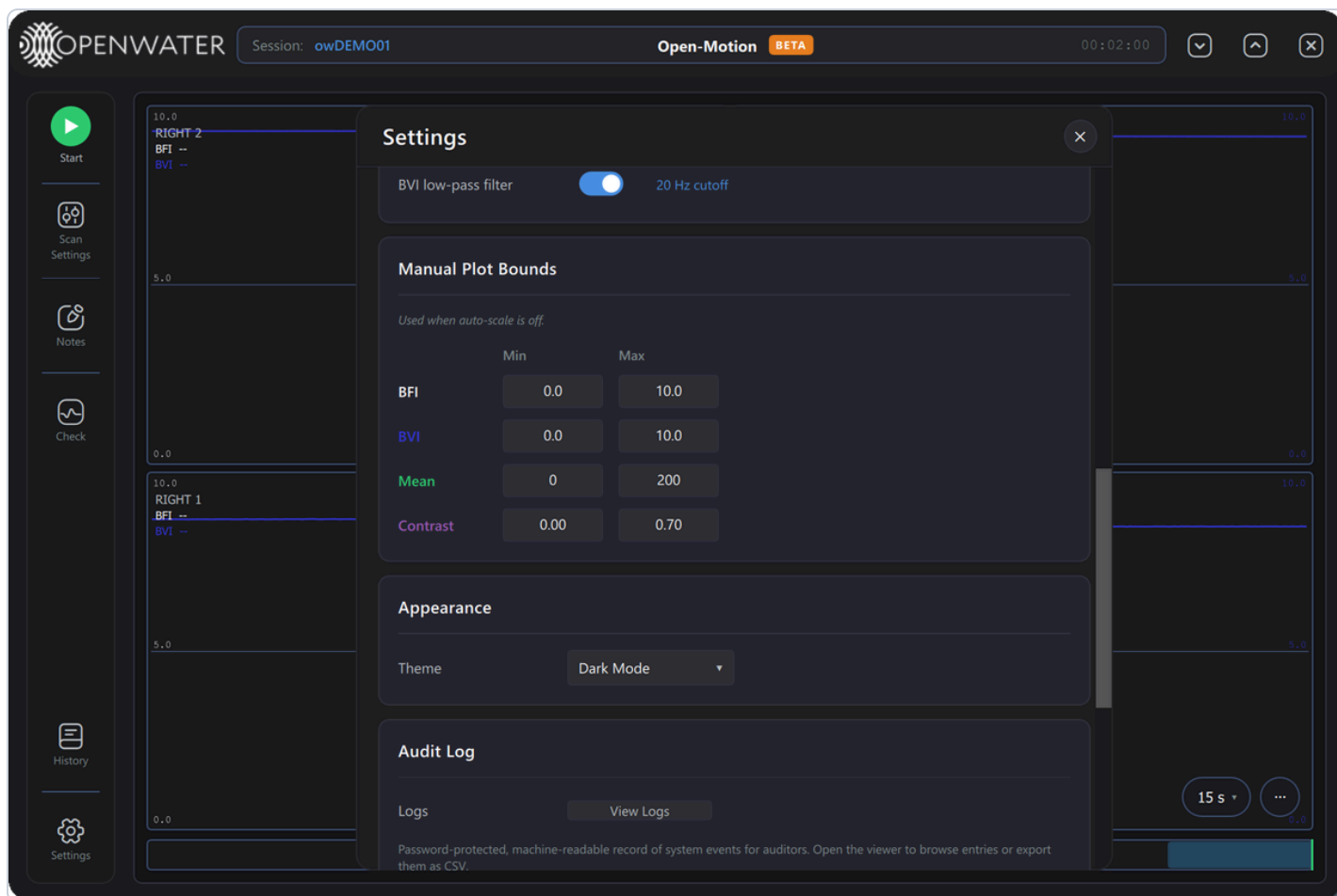

Card	Contents
Manual Plot Bounds	Fixed vertical ranges used when autoscale is off: Min and Max for BFI and BVI (0–10), Mean (0–1024) and Contrast (0–1). Min always stays below Max; out-of-range entries are corrected automatically.
Appearance	Theme: Dark Mode (default), Light Mode , or Liquid Glass . The interface re-themes immediately.
Audit Log	View Logs opens the audit log — a machine-readable record of system events for auditors. It asks for a password issued by Openwater; the viewer filters by event type, date range and text and exports CSV.

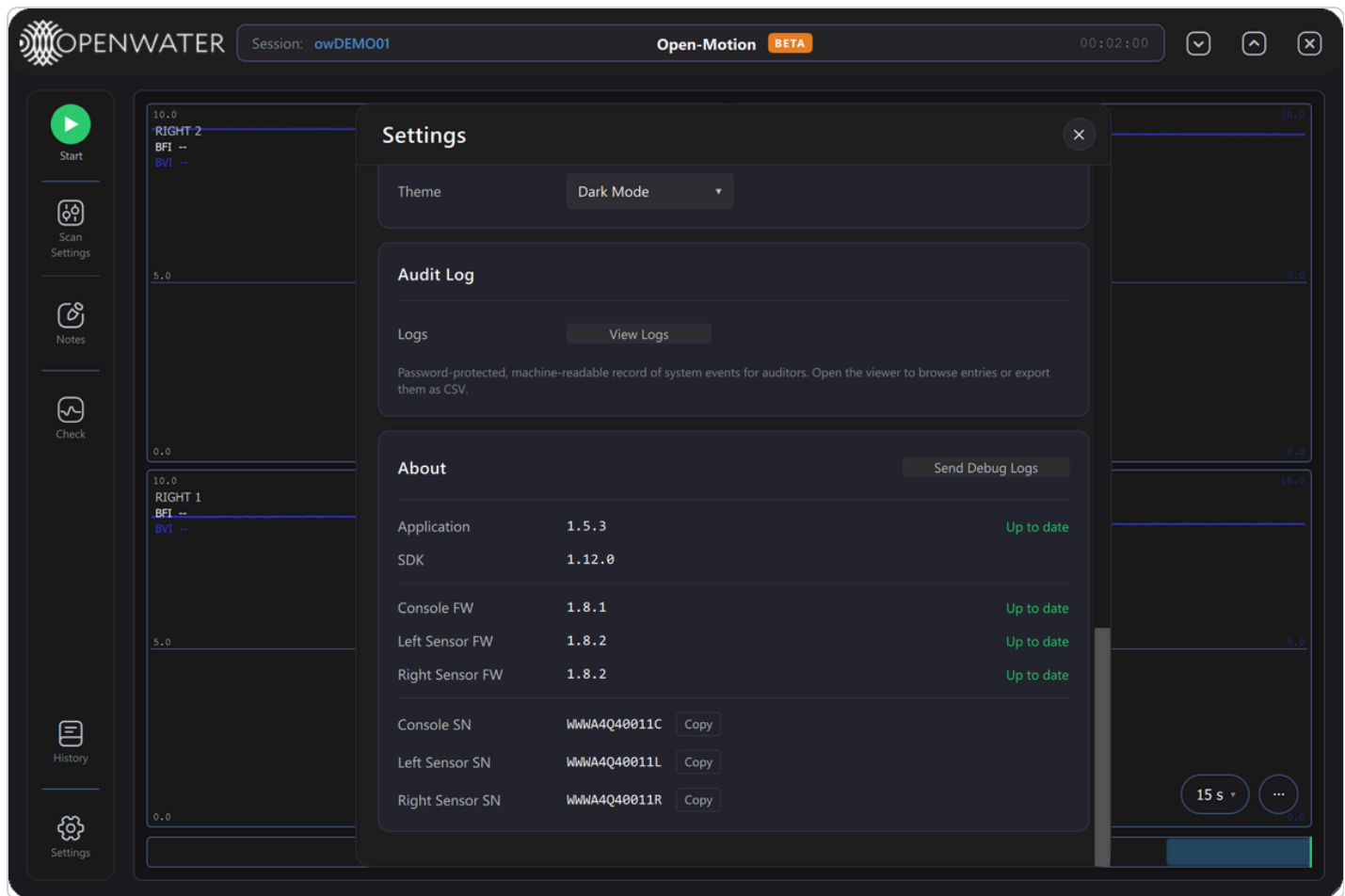

The **About** card lists:

Row	Meaning
Application	The application version, with <i>Up to date</i> , <i>Check failed</i> (no internet access or GitHub unreachable), or an Update button (§2.2).
SDK	The version of the Open-Motion software library.
Console FW, Left Sensor FW, Right Sensor FW	Firmware version of each connected device (or <i>Not connected</i>), with <i>Up to date</i> or an Update button (§2.3).
Console SN, Left Sensor SN, Right Sensor SN	Hardware serial numbers. Copy puts a serial number on the clipboard — include them when contacting support.
Send Debug Logs	Packages recent application logs for support (§13).

11. Data storage

Each Windows account keeps its own data in a private folder:

```
C:\Users\<Windows user>\AppData\Local\Openwater
```

Location	Contents
logs\	One application log file per launch (open-motion-<date>_<time>.log).
data\scans.db	The scan database: all scans, their notes, and your display preferences.
data\<full label>_left_mask<xx>_raw.csv , ..._right_mask<xx>_raw.csv	Raw camera histograms, one file per sensor, only when <i>Save raw</i> CSV is on. These are large — several hundred megabytes per camera per hour.
data\debug-bundles\	Zip files created by <i>Send Debug Logs</i> .
data\updates\	Installers downloaded by the in-app updater.

- Scan data is kept automatically in the database. Processed CSV files are only created when you export from Scan History (§9.1).
- In the Research build the database is not encrypted. Protect the computer and any exported files according to your study's data-handling rules.
- Scans recorded by version 1.5.2 or earlier remain in `C:\ProgramData\Openwater\data` and are not shown by this version (§2.1).
- On macOS, the data folder is `~/Library/Application Support/Openwater` .

12. Messages and errors

12.1 Messages (toasts)

Short messages appear in the bottom-right corner, color-coded: **green** = success, **yellow** = warning, **red** = error, **blue** = information. Most dismiss themselves after a few seconds (hovering over one pauses its countdown); those with an  can be closed by hand. Up to five are shown at a time.

Message	Meaning / action
"Note saved."	Session notes were stored.
"Console not detected...", "Sensor not detected...", "System not found..."	A device was not found at startup (§2.4). Check cables and power.
"Left sensor disconnected" (or <i>Right sensor, Console</i>)	A device was unplugged or powered off while no scan was using it. Reconnect it.
"Could not start scan — the previous step is still finishing. Please press Start again."	The system was still busy. Wait a few seconds and press Start again.
"Scan failed: ..."	The scan could not start or run; the reason follows.
"Scan duration cannot be 0 seconds — reset to 1 minute."	See §4.3.
"Camera ... connection lost ...", "... reconnected ...", "... temperature ..."	Individual camera events during a scan (§6).
"Scan ended unexpectedly ..."	See §6.
"Laser safety warning detected. Please restart your console. If this error persists, please contact support." (red, stays on screen)	The console's hardware laser-safety monitor reported a fault. Power-cycle the console. Contact support if it recurs.
"Exported to ...", "Exported N scan(s) to ..."	A CSV export finished.
"Interrupted scans can't be exported."	The selected scan was interrupted and has no exportable data.
"Console SN copied to clipboard." (or <i>Left / Right Sensor SN</i>)	A serial number was copied from Settings → About.
"Preparing debug logs...", "Debug logs saved to ..."	See §13.

12.2 Closing while busy

Pressing  during a scan or check shows "... in progress. Click X again to cancel and exit." Press  again within 5 seconds only if you intend to abandon the work.

12.3 Error dialog

Serious conditions raise a blocking **ERROR** dialog. It shows the error code (for example `E-202`), a title, an explanation, and a suggested action (→). While it is shown, the console light blinks blue.

Control	What it does
► Details	Expands technical detail, when available.
N more	Several errors are queued; they are shown one after another.
Copy details	Copies the error report to the clipboard.
Contact Support	Opens your email program with a pre-filled message to Openwater support, copies the report to the clipboard, and opens the folder containing the application log with the file selected — attach it to the email before sending.
Dismiss	Closes the dialog (or shows the next queued error). Clicking outside the dialog or pressing <code>Esc</code> does the same.

Quote the error code when contacting support.

12.4 Error code reference

Code	Title	Meaning	What to do
<code>E-101</code>	Sensor self-check failed	One or more of a sensor module's internal devices did not respond to its power-on self-check. The system cannot scan reliably.	Power-cycle the sensor module and reconnect. If it persists, the module needs service — contact support.
<code>E-102</code>	Sensor self-check unreadable	The sensor module connected but did not report its self-check result.	Power-cycle and reconnect. If it persists, update the sensor firmware or contact support.
<code>E-103</code>	Console initialization failed	The console's laser-power configuration could not be applied.	Power-cycle the console and reconnect. If it persists, contact support.
<code>E-104</code>	Console not detected	Startup check: no console found (yellow message, not a dialog).	Check the console USB cable and power.
<code>E-105</code>	Camera power-on failed	The sensor module could not power on its cameras.	Power-cycle the sensor module and reconnect. If it persists, the module may need service.
<code>E-106</code>	Sensor not detected	Startup check: no sensor module found (yellow message, not a dialog).	Check the sensor USB cable and power.

Code	Title	Meaning	What to do
E-201	Laser safety monitor unresponsive	Laser safety could not be confirmed, so the laser was shut off as a precaution.	Power-cycle the system and reconnect. Do not scan until this clears. Contact support if it persists.
E-202	Laser safety trip	The laser-safety monitor tripped during a scan. The laser was shut off and the scan stops about 5 seconds later.	Remove any obstruction, let the system settle, and start a new scan. If it trips repeatedly, stop and contact support.
E-301	Scan aborted before start	A check before the laser fired failed, so the scan did not start.	Resolve the reported problem (connection, safety or configuration) and start again.
E-302	Scan could not start	The system refused a new scan, usually because the previous one is still finishing.	Wait a few seconds and try again. Reconnect if it stays busy.
E-303	Camera data lost during scan	Every camera stopped delivering data for about 3 seconds, so the scan was stopped. Data captured before the loss was saved.	Check the sensor cables and power, reconnect, and start a new scan.
E-304	Device disconnected during scan	The console or a sensor module in use was disconnected during the scan, so the scan was stopped. Data captured before the disconnection may not have been saved.	Check the USB cables and power, reconnect, and start a new scan.

13. Getting help

1. Open **Settings** → **About** and press **Send Debug Logs**. The application packages the last 48 hours of application logs into a zip file under `data\debug-bundles\`, opens the folder with the file selected, and shows the file's location. The package contains application logs only — no scan data.
2. Email the zip file to Openwater support at **support@openwater.health** with a description of the problem, the error code if there was one, and the serial numbers from Settings → About.
3. From an error dialog, **Contact Support** prepares the email for you (§12.3).

14. Known issues in 1.5.3

Issue	Workaround
After changing the User Label between scans, the data file name of the next scan can still carry the previous label; it is correct from the scan after that.	Check the file names after changing the label.
The "Debug logs saved..." message names <code>support@openwater.cc</code> .	Send debug logs to the support address in §13.
Plot data loads slowly when zooming out to the beginning of a long scan.	Allow a moment for the plot to fill in.
A scan just under a whole minute can show a duration such as <code>0:60</code> in Scan History.	Read it as <code>1:00</code> . The duration line in the scan's notes is correct.
When the window is maximized, the restore button shows an incorrect icon.	The button still restores the window.